



STERISIL®
Safe, Simple, Sterisil

○ SHOCK
○ TEST
✓ MAINTAIN

Sterisil® Cartridge

THE MOST POWERFUL DIRECT FEED WATERLINE
TREATMENT CARTRIDGES ON THE MARKET



*Quick disconnect kit, shown installed
(included)*



Valved Cartridge setup kit (sold separately)

KEY FEATURES

- Delivers purified and treated dental water direct to dental units
- EPA validated treatment claim: ≤ 10 CFU/ml
- Shock treatment included
- Prevents hard water deposits
- Valved and inline configurations available
- Safe for patients, staff, and equipment
- Does not adversely affect composite bond strength
- Amalgam separator safe
- Multiple sizes available

The Sterisil® Cartridge delivers purified and treated dental water to chairs without a bottle. Utilizing the same ion exchange technology found in the Sterisil® Straw V2, Cartridges include shock and maintenance treatments and last up to 1 year*.

Visit us online at Sterisil.com

*Cartridge lifespans vary with cartridge size and water quality.
Water quality is assessed in parts per million (PPM) total dissolved solids (TDS).





STERISIL®
Safe, Simple, Sterisil

CORE FEATURES THAT DELIVER COMPLIANT DENTAL WATER

SHOCK TREATMENT

The initial shock treatment releases automatically when water is introduced.

CONTINUOUS MAINTENANCE

Once installed on the dental unit, the cartridge is continuously dispensing residual silver, keeping bacteria in check.

SIMPLE PURIFICATION

Variable tap water quality is no longer a concern. Filter media inside the cartridge purifies water before it's treated with silver.

VALVED OR INLINE


Choose between valved or inline models in multiple cartridge sizes for the perfect fit regardless of your space constraints.

USE THE LARGEST CARTRIDGE THAT WILL FIT WITHIN YOUR SPACE CONSTRAINTS

STERISIL PART #	DESCRIPTION
Inline Cartridges (setup kit included)	
Ci-3i	Length 6-1/16" X diameter 2"
Ci-7i	Length 8-1/16 X diameter 2.5"
Ci-9i	Length 10-1/8" X diameter 2.5"
Valved Cartridges (setup kit NOT included)	
CV-10	Length 10-3/8" X diameter 2.5"
CV-20	Length 14-3/4" X diameter 2.5"
CV-SK	Fittings, valve, head and mounting bracket for valved cartridges




Proudly made
in the USA,

 Sterisil, Inc.
835 S. Hwy 105 Ste. D
Palmer Lake, CO USA 80133
+1 719 622 7200

EPA Reg. No. 83315-2
EPA Est. No. 83315-CO-001



1639

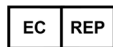

EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands



For more information contact your **Sterisil® sales representative or local dealer.**

719 622 7200 or online at Sterisil.com

GP4000.D

