



SHOCK TEST MAINTAIN

FASTCheck15™ In-Office Water Test Strip

The fastest, most affordable way to monitor water quality for harmful bacteria




FAST

- Quick, 15-minute in-office water bacteria test
- Detects microbials exceeding 500CFU/mL
- Know your water is safe in 15 minutes vs. 3 days for paddle test results
- FASTCheck15™ removes the risk of bacteria growing to harmful levels in between water sampling and results

SIMPLE TO USE. SIMPLE TO READ.

- 4 simple steps to reduce user error
- No incubation/neutralizer needed
- No counting red dots or color-matching
- Easy-to-read test strip that provides clear results: 1 line = Pass, 2 lines = Fail

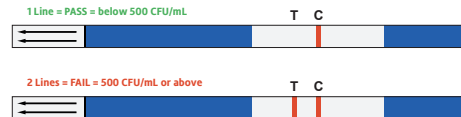

SCAN FOR SIMPLE TEST INSTRUCTIONS AND VIDEO
[Bit.ly/fastcheck15-ifu](https://bit.ly/fastcheck15-ifu)

LOWEST COST WITH HIGH ACCURACY

- A fraction of the cost of other waterline tests
- 95% accurate results



My Solutions Center™

SCAN TO DOCUMENT TEST RESULTS
Sterisil recommends following the CDC, OSAP, ADA, and OSHA guidelines to document water monitoring results. Find test logs and support at My Solutions Center.
[Bit.ly/mysolutionscenter](https://bit.ly/mysolutionscenter)

Solmetex System of Solutions

Source Water Quality



Dental Unit Waterline Management



Isolation and Aerosol Control



Evacuation




Amalgam Capture and Recycling



Waste Management



Every Drop Matters™

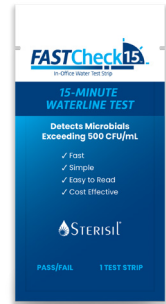

The Fastest and Most Convenient

Stay ahead of the biofilm problem and eliminate having an empty chair, even if you fail a test. Easy-to-read results in just 15 minutes makes it possible to test after your last patient on any day, and be ready for your next patient the following morning. FASTCheck15™ gives you control over your waterline safety in minutes, not days.

Be proactive about your waterlines instead of reactive. Passing more tests mean your staff can spend their time on what's most important to you... your patients.

Sterisil® SafeWater Testing Protocol

Stay ahead of biofilm by testing and shocking monthly.



Month 1

Test in-office with FASTCheck15™



Shock with Citrisil™



Month 2

Test in-office with FASTCheck15™



Shock with Citrisil™



Month 3

Mail-in test with Sterisil® R2A™ Waterline test



Shock with Citrisil™



1-2-3 REPEAT for proactive waterline health



Scan to document results at:



my-solutions-center.solmetex.com

Failed test? It happens. Here's what to do:

Retest 5 days after shocking* If retest results fail, contact our experienced Customer Support specialists at **800-216-5505** or chat with a Water Specialist online at **Solmetex.com**

*Experts recommend waiting at least 5 days before retesting.

Testing Protocol FAQs

solmetex.com/failedtestprotocol



*CDC, OSAP, FDA, and ADA recommend documenting water monitoring results.

STERISIL PART #	SHELF LIFE	DESCRIPTION
FC15-6	2 YEARS	PACK OF 6 FASTCHECK15™ TEST STRIPS
FC15-1	2 YEARS	1 SAMPLE FASTCHECK15™ TEST STRIP (Not available for retail.)

Order/Info:

(800) 216-5505
Solmetex.com

SHOP NOW!
[Bit.ly/fastcheck15-facts](https://bit.ly/fastcheck15-facts)



Solmetex™
Dentistry's Complete
Provider of Water Solutions

Solmetex
Sterisil
DryShield